



Product Sheet

# *Methanothermobacter thermautotrophicus* (ATCC® 29183™)

Please read this **FIRST**



## Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

## Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Methanothermobacter thermautotrophicus* (ATCC® 29183™)

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor

## Description

Designation: YT

Deposited Name: *Methanothermobacter thermoautotrophicum* Zeikus and Wolfe

## Propagation

### Medium

ATCC® Medium 1045: Methanobacteria medium

### Growth Conditions

Temperature: 60.0°C

Atmosphere: Under a gas mixture of 80% H<sub>2</sub>, 20% CO<sub>2</sub>

### Propagation Procedure

1. Sterilize the top of the Balch tube by spraying it with 70% ethanol and then flame the top.
2. If needed exchange the gas in the test tube for 80% H<sub>2</sub> 20% CO<sub>2</sub>.
4. When the Balch tube is ready to inoculate, thaw the frozen vial at room temperature under a gentle stream of oxygen-free gas.
5. For inoculation, use an anaerobic (see c below) 1.0 ml syringe tipped with 22-gauge needle, withdraw the cell suspension from the vial and transfer it to the broth. Plate 0.1 ml of the inoculated culture onto a non-selective medium and incubate aerobically at 37°C. Use 0.1 ml of the inoculated culture to inoculate a nonselective aerobic broth and an additional tube of #1045 broth. Incubate the non-selective aerobic broth tubes at 37°C. Incubate the #1045 broth culture at 60°C.
6. Growth should be detected in the #1045 broth within 24 to 48 hours. There should be no growth detected on the aerobic plate or in the aerobic broth.

### ANAEROBIC CONDITIONS:

- a. Resazurin is a commonly used redox indicator that is pink when the redox potential is above 50 mv., and colorless when the redox potential is below 110 mv. i.e. highly reducing. Most strict anaerobes require this low redox potential for optimum growth.
- b. To obtain a fully reduced medium, it is necessary that the medium be anoxic and that a reducing agent be added. Common reducing agents are sodium sulfide, cysteine, dithiothreitol, and titanium citrate.
- c. Syringes can be made anaerobic by one of two methods. 1. Displace the dead space in the syringe with a sterile

## Notes

At 1000x magnification cells are long rods that occur singly, pairs, and in clumps.

## References

References and other information relating to this product are available online at [www.atcc.org](http://www.atcc.org).

## Biosafety Level: 1

Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes for Health.

## ATCC Warranty

The viability of ATCC® products is warranted for 30 days from the date of shipment, and is valid only if the product is stored and cultured according to the information included on this product information sheet. ATCC lists the media formulation that has been found to be effective for this strain. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this strain. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

## Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans.



Product Sheet

***Methanothermobacter  
thermautotrophicus***  
**(ATCC® 29183™)**

While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.


This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, ATCC is not liable for damages arising from the misidentification or misrepresentation of cultures.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this product. The MTA is also available on our Web site at [www.atcc.org](http://www.atcc.org)

Additional information on this culture is available on the ATCC web site at [www.atcc.org](http://www.atcc.org).  
© ATCC 2013. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [02/18]

---

**Please read this FIRST**

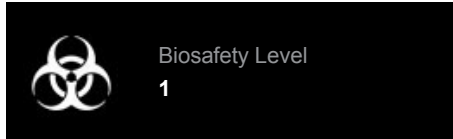

---

**Intended Use**

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

---

**Citation of Strain**

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Methanothermobacter thermautotrophicus* (ATCC® 29183™)

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor