




Product Sheet

# *Paramecium primaurelia* (ATCC® 30756™)

Please read this FIRST

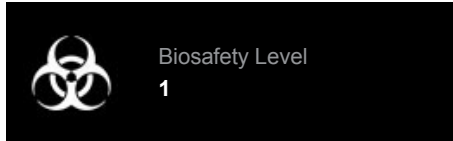

## Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

## Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Paramecium primaurelia* (ATCC® 30756™)

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor

## Description

**Strain Designation:** stock 183  
**Deposited Name:** *Paramecium primaurelia* Sonneborn  
**Depositor:** TM Sonneborn  
**Isolation:**  
Tokyo, Japan, pre-1951

## Propagation

### Growth Conditions

**Max Temperature:** 31.0°C

**Min Temperature:** 13.0°C

**Protocol:** ATCCNO: 30300 SPEC: This strain is shipped as a growing test tube culture. Upon arrival, remove test tube from sealed plastic envelope, remove plastic seal from cap, and loosen the cap one half turn. Add 1.0 ml of ATCC medium 802 bacterized with *Enterobacter aerogenes* ATCC 13048 twice weekly. When the tube is filled to within one inch of the top, decant leaving 5.0 ml in the original tube. Subcultures are established by transferring 0.5 ml of a growing culture to 5.0 ml of bacterized ATCC medium 802 in a 20 x 120 mm test tube.

### Medium

ATCC® Medium 802: Sonneborn's *Paramecium* medium

### Instructions for Complete Medium

ATCC Medium 802 inoculated with *Enterobacter aerogenes* (ATCC® 13048)

### Culture Maintenance

Subculture every two months to a fresh tube of bacterized medium in the following manner:

1. Transfer 0.5 ml from a growing culture to 5.0 ml of ATCC medium 802 bacterized with *Enterobacter aerogenes* (ATCC® 13048).
2. Add 1.0 ml of ATCC medium 802 bacterized with *Enterobacter aerogenes* ATCC 13048 twice weekly. When the tube is filled to within one inch of the top, aspirate from the bottom of the tube and reduce the volume to 5.0 ml.
3. Incubate upright at 25°C with the caps on loosely.

## Cryopreservation

### Cryoprotective Solution

DMSO	1.5 ml
Fresh growth medium w/o bacteria	7.5 ml
MgCl <sub>2</sub> (0.5 mM)	0.5 ml
CaCl <sub>2</sub> (0.5 mM)	0.5 ml


1. Mix the components in the order listed. Before adding the MgCl<sub>2</sub> and the CaCl<sub>2</sub> allow the solution to return to room temperature. When the medium is added to the DMSO the solution will warm up due to chemical heat.
2. Harvest cells from a culture that is at or near peak density by filtration and centrifugation at 200 x g for 1 min.
3. Adjust the concentration of cells to 2 x 10<sup>5</sup>/ml in fresh medium.
4. Mix the cell preparation and the cryoprotective solution in equal portions.
5. Dispense in 0.5 ml aliquots into 1.0 - 2.0 ml sterile plastic screw-capped cryules (special plastic vials for cryopreservation).
6. Place vials in a controlled rate freezing unit. From room temperature cool at -1°C/min to -40°C. If freezing unit can compensate for the heat of fusion, maintain rate at -1 C/min through heat of fusion. At -40°C plunge ampules into liquid nitrogen.
7. Ampules are stored in either the vapor or liquid phase of a nitrogen refrigerator.
8. To establish a culture from the frozen state add 1.0 ml ATCC medium 802 to the frozen ampule and place it in a 35°C water bath. Immerse the vial to a level just above the surface of the frozen material. Do not agitate the vial.
9. Immediately after thawing, do not leave in water bath, aseptically remove the contents of the ampule and inoculate onto the surface of an ATCC medium 919 (non-nutrient agar) plate containing an overlay of 15.0 ml of bacterized ATCC medium 802.
10. Incubate at 25°C.



Product Sheet

## *Paramecium primaurelia* (ATCC® 30756™)

Please read this **FIRST**

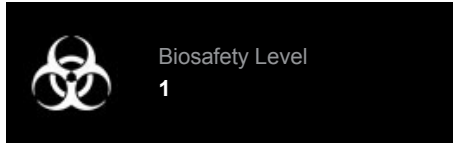

### Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

### Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Paramecium primaurelia* (ATCC® 30756™)

11. Once the culture is established, transfer 0.5 ml to 5.0 ml of bacterized ATCC medium 802.
12. Follow the protocol for maintenance of culture.

### References

References and other information relating to this product are available online at [www.atcc.org](http://www.atcc.org).

### Biosafety Level: 1

Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes for Health.

### ATCC Warranty

The viability of ATCC® products is warranted for 30 days from the date of shipment, and is valid only if the product is stored and cultured according to the information included on this product information sheet. ATCC lists the media formulation that has been found to be effective for this strain. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this strain. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

### Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans.

While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, ATCC is not liable for damages arising from the misidentification or misrepresentation of cultures.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this product. The MTA is also available on our Web site at [www.atcc.org](http://www.atcc.org)

Additional information on this culture is available on the ATCC web site at [www.atcc.org](http://www.atcc.org).

© ATCC 2013. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [02/06]

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor