



Product Sheet

## *Diplonema ambulator* (ATCC® 50223™)

Please read this **FIRST**



### Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

### Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Diplonema ambulator* (ATCC® 50223™)

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor

### Description

**Strain Designation:** X3Smith

**Deposited Name:** *Diplonema ambulator* Larsen and Patterson

**Depositor:** TA Nerad

**Isolation:**

artificial Chesapeake Bay Marine Ecosystem, Smithsonian Institution, Washington, DC, 1987

### Propagation

#### Growth Conditions

**Temperature:** 25.0°C

Duration: axenic

#### Medium

ATCC® Medium 1728: Enriched Isonema medium

#### Culture Maintenance

1. Vigorously agitate a culture at or near peak density and aseptically transfer 0.1-0.2 ml to 10 ml of fresh ATCC medium 1728 in a T-25 tissue culture flask.
2. Incubate in a 25°C incubator with the caps screwed on tightly.

### Cryopreservation

1. Harvest the culture by agitating the contents of each flask. Transfer the cell suspensions to 15 ml plastic centrifuge tubes.
2. Spin the cell suspensions at approximately 500 x g for 5 min.
3. Pool the cell pellets and adjust the concentration to 2.0 - 4.0 x 10<sup>7</sup> cells/ml with a fresh ATCC medium 1728.  
\*If the concentration is too low centrifuge at 500 x g for 5 min and resuspend in the volume of ATCC medium 1728 required to yield the desired concentration.
4. Mix the cell preparation and 20% (v/v) DMSO in equal portions. The final concentration will be 1.0 - 2.0 x 10<sup>7</sup> cells/ml and 10% DMSO. The time from the mixing of the cell preparation and cryoprotective solution to the start of the freezing process should be no less than 15 min. and no more than 30 min.
5. Dispense in 0.5 ml aliquots to 1.0-2.0 ml sterile plastic screw-capped cryules (special plastic vials for cryopreservation).
6. Place the vials in a controlled rate freezing unit. From room temperature cool at -1°C/min to -40°C. If the freezing unit can compensate for the heat of fusion, maintain rate at -1°C/min through the heat of fusion. At -40°C plunge into liquid nitrogen. Alternatively, place the vials in a Nalgene 1°C freezing apparatus. Place the apparatus at -80°C for 2.5 to 3 hours and then plunge ampules into liquid nitrogen. (The cooling rate in this apparatus is approximately -1°C/min.)
7. Store in either the vapor or liquid phase of a nitrogen refrigerator.
8. To thaw a frozen ampule, place it in a 35°C water bath such that the lip of the ampule remains above the water line. Thawing time is approximately 2 to 3 minutes. Do not agitate the ampule. Do not leave ampule in water bath after thawed.
9. Immediately after thawing, aseptically transfer contents to a T-25 tissue culture flask containing 10 ml ATCC medium 1728.
10. Incubate in a 25°C CO<sub>2</sub> incubator with the caps screwed on tightly.

### References

References and other information relating to this product are available online at [www.atcc.org](http://www.atcc.org).

### Biosafety Level: 1

Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes for Health.




Product Sheet

## *Diplonema ambulator* (ATCC® 50223™)

---

### Please read this FIRST

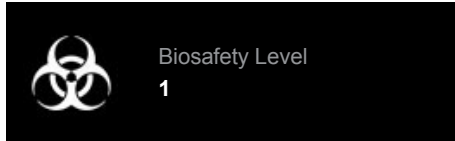

---

### Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

---

### Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Diplonema ambulator* (ATCC® 50223™)

---

### ATCC Warranty

The viability of ATCC® products is warranted for 30 days from the date of shipment, and is valid only if the product is stored and cultured according to the information included on this product information sheet. ATCC lists the media formulation that has been found to be effective for this strain. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this strain. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

---

### Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans.

While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, ATCC is not liable for damages arising from the misidentification or misrepresentation of cultures.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this product. The MTA is also available on our Web site at [www.atcc.org](http://www.atcc.org)

Additional information on this culture is available on the ATCC web site at [www.atcc.org](http://www.atcc.org).

© ATCC 2013. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [02/06]

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor