



Product Sheet

Dialister pneumosintes (ATCC® 51894™)

Please read this **FIRST**



Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Dialister pneumosintes* (ATCC® 51894™)

Description

Designation: N20A-27

Deposited Name: *Dialister pneumosintes* (Olitsky and Gates) Moore and Moore

Propagation

Medium

ATCC® Medium 1490: Modified chopped meat medium

Growth Conditions

Temperature: 37.0°C

Atmosphere: Anaerobic

Propagation Procedure

1. Open vial according to enclosed instructions.
2. Under anaerobic conditions, withdraw 0.5 ml of recommended broth from a single test tube (5 to 6 ml) and rehydrate the entire vial contents.
3. Aseptically transfer this aliquot back into the broth tube. A slant and a pre-reduced blood plate may also be inoculated with 0.1 ml each of the cell suspension. An aerobic blood plate may also be streaked to check for purity.
4. Incubate tubes and one plate under anaerobic conditions at 37°C. Incubate one blood plate aerobically at 37°C.
5. After 4 to 5 days, colonies are pinpoint and transparent. Growth in broth shows only slight turbidity and a smooth sediment.

ANAEROBIC CONDITIONS:

Anaerobic conditions for transfer may be obtained by either of the following:

- Use of an anaerobic gas chamber, or
- Placement of test tubes under a gassing cannula system hooked to anaerobic gas.
Anaerobic conditions for incubation may be obtained by any of the following:
 - Loose screw caps on test tubes in anaerobic chamber,
 - Loose screw caps on test tubes in an activated anaerobic gas pack jar, or
 - Use of sterile butyl rubber stoppers on test tubes so that an anaerobic gas headspace is retained.

References

References and other information relating to this product are available online at www.atcc.org.

Biosafety Level: 1

Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes for Health.

ATCC Warranty

The viability of ATCC® products is warranted for 30 days from the date of shipment, and is valid only if the product is stored and cultured according to the information included on this product information sheet. ATCC lists the media formulation that has been found to be effective for this strain. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this strain. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans.

While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, ATCC is not liable for

American Type Culture Collection
PO Box 1549
Manassas, VA 20108 USA
www.atcc.org

800.638.6597 or 703.365.2700
Fax: 703.365.2750
Email: Tech@atcc.org


Or contact your local distributor



Product Sheet

Dialister pneumosintes **(ATCC® 51894™)**

Please read this FIRST

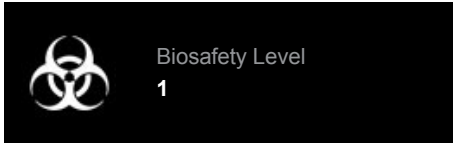

Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Dialister pneumosintes* (ATCC® 51894™)

damages arising from the misidentification or misrepresentation of cultures.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this product. The MTA is also available on our Web site at www.atcc.org

Additional information on this culture is available on the ATCC web site at www.atcc.org.

© ATCC 2013. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [02/02]

American Type Culture Collection
PO Box 1549
Manassas, VA 20108 USA
www.atcc.org

800.638.6597 or 703.365.2700
Fax: 703.365.2750
Email: Tech@atcc.org

Or contact your local distributor