




Product Sheet

Chlamydomonas reinhardtii (ATCC® PRA-141™)

Please read this **FIRST**

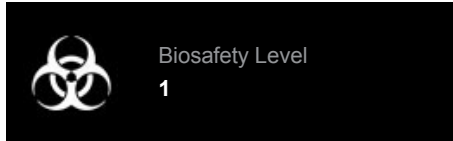

Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Chlamydomonas reinhardtii* (ATCC® PRA-141™)

American Type Culture Collection
PO Box 1549
Manassas, VA 20108 USA
www.atcc.org

800.638.6597 or 703.365.2700
Fax: 703.365.2750
Email: Tech@atcc.org

Or contact your local distributor

Description

Strain Designation: CC-1690 wt mt+ Sager 21 gr
Deposited Name: *Chlamydomonas reinhardtii* Dangeard
Depositor: EH Harris, GM Smith
Isolation:
soil
Amherst Massachusetts, United States
Isolation date: December, 1945

Propagation

Growth Conditions

Duration: axenic
Protocol: Streak onto agar slants
Interval: 3 to 6 months on agar; weekly in liquid

Instructions for Complete Medium

ATCC Medium 5

Protocols

Frozen ampules packed in dry ice should either be thawed immediately or stored in liquid nitrogen. If liquid nitrogen storage facilities are not available, frozen ampules may be stored at or below -70°C for approximately one week. **Do not under any circumstance store frozen ampules at refrigerator freezer temperatures (generally -20°C).** Storage of frozen material at this temperature will result in the death of the culture.

1. To thaw a frozen ampule, place in a 35°C water bath, until thawed (2-3 min). Immerse the ampule just sufficient to cover the frozen material. Do not agitate the ampule.
2. Immediately after thawing, aseptically transfer the entire contents to a single 16 x 125 mm screw-capped test tube containing 5 ml of ATCC Medium 5 broth. Incubate the tube on a 15° horizontal slant with the cap screwed on loosely (loosened one half turn) at 25°C under a 14 hour light (~50 mEinsteins/m²/s irradiance)/10 hour dark cycle. Alternatively, add the entire thawed contents to the surface of a 20 x 100 mm Petri plate containing 20 ml of ATCC medium 5 agar. Wrap the plate culture with parafilm and incubate upright under the same light/dark cycle as specified for a test tube culture.

Culture Maintenance

1. Screw the cap on tightly and vigorously agitate the culture.
2. Aseptically transfer a 0.1ml aliquot to 5 ml of fresh medium in a 16 x 125 mm screw-capped test tube.
3. Screw caps on loosely (loosened one-half turn) and incubate on a 15° horizontal slant at 25°C under a 14 hour light (~50 mEinsteins/m²/s irradiance)/10 hour dark cycle.
4. Subculture every 14-21 days.

Cryopreservation


1. Harvest cells from a culture that is at or near peak density by centrifugation at 800 x g for 5 min.
2. Adjust the concentration of cells to 2 x 10⁶ - 2 x 10⁷/ml in fresh medium.
3. While cells are centrifuging prepare a 10% (v/v) solution of sterile methanol in fresh medium.
4. Mix the cell preparation and the 10% methanol in equal portions. Thus, the final concentration will be 10⁶ - 10⁷ cells/ml and 5% (v/v) Methanol. The time from the mixing of the cell preparation and methanol stock solution to the beginning of the freezing process should be no less than 5 min and no greater than 15 min.
5. Dispense in 0.5 ml aliquots into 1.0 - 2.0 ml sterile plastic screw-capped cryovials (special plastic vials for cryopreservation).
6. Place the vials in a controlled rate freezing unit. From room temperature cool at -1°C/min to -40°C. If the freezing unit can compensate for the heat of fusion, maintain rate at -1°C/min through the heat of fusion. At -40°C plunge into liquid nitrogen. Alternatively, place the vials in a Nalgene 1°C freezing apparatus. Place the apparatus at -80°C for 1.5 to 2 hours and then plunge ampules into liquid nitrogen. (The cooling rate in this apparatus is approximately -1°C/min.)
7. The frozen preparations should be stored in either the vapor or liquid phase of a nitrogen refrigerator.



Product Sheet

Chlamydomonas reinhardtii (ATCC® PRA-141™)

Please read this FIRST

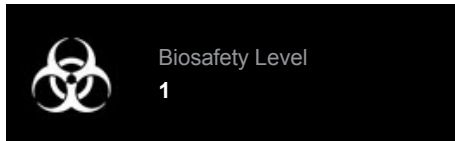

Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Chlamydomonas reinhardtii* (ATCC® PRA-141™)

- Frozen preparations stored below -130°C are stable indefinitely. Those stored at temperatures above -130°C are progressively less stable as the storage temperature is elevated. Vials should not be stored above -55°C.
- To establish a culture from the frozen state place an ampule in a water bath set at 35°C. Immerse the vial just to a level just above the surface of the frozen material. Do not agitate the vial.
 - Immediately after thawing, do not leave in the water bath, aseptically remove the contents of the ampule and add to a centrifuge tube containing 5 ml of ATCC medium 5 without agar. Centrifuge at 300 x g for 5 min.
 - Remove most of the supernatant (=methanol, which can inhibit growth) and then resuspend the pellet. Transfer the culture to a 16 x 125 mm screw-capped test tube containing 5 ml of ATCC medium 5 broth or to the surface of an ATCC medium 5 agar plate (20 x 100 mm Petri plate containing 20 ml of ATCC medium 5 agar).
 - Incubate the culture at 50-100 µEinsteins/m²/s irradiance at 25°C. Maintain under a 14/10h light-dark photoperiod.



References

References and other information relating to this product are available online at www.atcc.org.



Biosafety Level: 1

Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes for Health.

ATCC Warranty

The viability of ATCC® products is warranted for 30 days from the date of shipment, and is valid only if the product is stored and cultured according to the information included on this product information sheet. ATCC lists the media formulation that has been found to be effective for this strain. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this strain. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans.

While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, ATCC is not liable for damages arising from the misidentification or misrepresentation of cultures.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this product. The MTA is also available on our Web site at www.atcc.org

Additional information on this culture is available on the ATCC web site at www.atcc.org.

© ATCC 2013. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [02/26]

American Type Culture Collection
PO Box 1549
Manassas, VA 20108 USA
www.atcc.org

800.638.6597 or 703.365.2700
Fax: 703.365.2750
Email: Tech@atcc.org

Or contact your local distributor