




## Nucleic Acid Product Sheet

# ***Salmonella enterica* subsp. *enterica* (ex Kauffmann and Edwards) Le Minor and Popoff serovar Typhimurium (ATCC® 13311D-5™)**

### Please read this FIRST

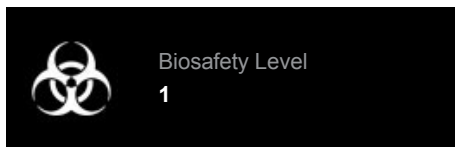

Biosafety Level  
1

### Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

### Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Salmonella enterica* subsp. *enterica* (ex Kauffmann and Edwards) Le Minor and Popoff serovar Typhimurium (ATCC® 13311D-5™)

### Nucleic Acid Information

Total DNA: At least 5 µg in 1X TE buffer.  
OD<sub>260</sub>/OD<sub>280</sub>: 1.6 to 2.0

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor

### Description

**Source:** *Salmonella enterica* subsp. *enterica* (ex Kauffmann and Edwards) Le Minor and Popoff serovar Typhimurium

**Designation:** Genomic DNA from *Salmonella enterica* subsp. *enterica* serovar Typhimurium strain NCTC 74 [ATCC® 13311™]

**Description:** Genomic DNA isolated from *Salmonella enterica* subsp. *enterica* serovar Typhimurium strain NCTC 74. This bacterial strain is also available as ATCC® Catalog No. 13311™.

**Note:** Genomic DNA isolated from bacteria is appropriate for PCR\* and other molecular biology applications. \*The polymerase chain reaction (PCR) process is covered by patents owned by Hoffmann-LaRoche Inc. Use of the PCR process requires a license.

### Batch-Specific Information

Refer to the Certificate of Analysis for batch-specific test results.

### Preparation Procedure

1. Thaw the vial at room temperature and immediately place on ice. Avoid exposing the DNA to repeated freeze-thaw cycles as it may result in degradation of the DNA.
2. Gently mix the sample to ensure an even distribution of material.
3. Briefly centrifuge the tube before opening to ensure all liquid is at the bottom.

### Quality Control Information

1. Bacterial genomic DNA is provided in a frozen form. Store at -20°C or lower (-80°C or vapor-phase liquid nitrogen (-196°C) is preferred for long-term storage, more than one month) Note: Do not store in freezers with a defrost cycle. This will expose the product to increased temperatures.
2. Total DNA by PicoGreen® measurement was found to be approximately 5 µg.
3. Purity: OD<sub>260</sub>/OD<sub>280</sub> ratio.
4. Integrity of DNA was determined by electrophoresis on a 1% agarose gel stained with SYBR Safe™, and was found to be of high molecular weight.
5. No RNA was detected by electrophoresis.
6. Functional activity was confirmed by PCR amplification 16S ribosomal RNA gene.
7. Identity confirmed by sequencing of 16S ribosomal RNA gene (first ~ 500 base pairs).

### Biosafety Level: 1

Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes for Health.

### ATCC Warranty

The viability of ATCC® products is warranted for 30 days from the date of shipment, and is valid only if the product is stored and cultured according to the information included on this product information sheet. ATCC lists the media formulation that has been found to be effective for this strain. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this strain. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

### Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans.

While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, ATCC is not liable for damages arising from the misidentification or misrepresentation of cultures.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this



product. The MTA is also available on our Web site at [www.atcc.org](http://www.atcc.org)

Additional information on this culture is available on the ATCC web site at [www.atcc.org](http://www.atcc.org).


© ATCC 2015. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [08/04]

## Nucleic Acid Product Sheet

# ***Salmonella enterica* subsp. *enterica* (ex Kauffmann and Edwards) Le Minor and Popoff serovar Typhimurium (ATCC® 13311D-5™)**

---

Please read this **FIRST**

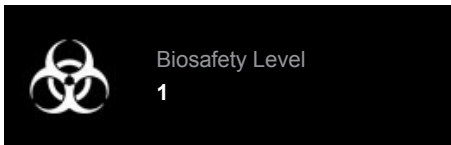

---

### Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

---

### Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Salmonella enterica* subsp. *enterica* (ex Kauffmann and Edwards) Le Minor and Popoff serovar Typhimurium (ATCC® 13311D-5™)

---

### Nucleic Acid Information

Total DNA: At least 5 µg in 1X TE buffer.

OD<sub>260</sub>/OD<sub>280</sub>: 1.6 to 2.0

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor