


Product Sheet

Toxoplasma gondii (ATCC[®] 50860™)

Please read this FIRST

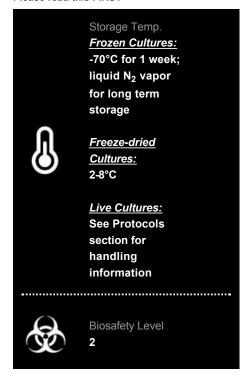

Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: $Toxoplasma\ gondii$ (ATCC® 50860TM)

American Type Culture Collection PO Box 1549 Manassas, VA 20108 USA www.atcc.org

800.638.6597 or 703.365.2700 Fax: 703.365.2750 Email: <u>Tech@atcc.org</u>

Or contact your local distributor

Description

Strain Designation: P110

Deposited Name: Toxoplasma gondii (Nicolle and Manceaux) Nicolle and Manceaux

Depositor: LD Sibley **Isolation:** Pig, USA **Genotype:** Haplogroup 3



Propagation

Growth Conditions
Temperature: 35°C to 37°C

Cell Line: ATCC® CRL-1634™ (human foreskin fibroblasts) (Contact ATCC Sales to order)

Instructions for Complete Medium

Dulbecco's Modified Eagle's Medium (DMEM) (ATCC® cat. 30-2002) with 4 mM L-glutamine adjusted to contain 1.5 g/L sodium bicarbonate, 4.5 g/L glucose, and 1.0 mM sodium pyruvate; supplemented with 3% fetal bovine serum (ATCC® cat. 30-2020). (Contact ATCC Sales to order)



Protocols

Cell Line Maintenance

- To establish a cell culture from the frozen state place an ampule in a water bath set at 35°C (2-3 min).
 Immerse the vial just sufficiently to cover the frozen material. Do not agitate the vial.
- Immediately after thawing, aseptically remove the contents of the ampule and inoculate into 10.0 mL of fresh ATCC® 30-2002 with 10% (v/v) Heat-Inactivated Fetal Bovine Serum (HIFBS)* in a T-25 tissue culture flask
- 3. Outgas the flask for 10 seconds with a 95% air, 5% CO₂ gas mixture.
- 4. Incubate in a 35-37°C CO₂ incubator with the cap screwed on tightly.
- 5. Change the medium 1-2 times per week.

*Fetal bovine serum is available from ATCC (catalog number 30-2020; contact ATCC Sales to order). Serum is heat-inactivated by exposure to 56°C for 30 minutes. This treatment will inactivate proteins of the complement pathway. Remove the serum from the refrigerator and aseptically distribute in 100 mL aliquots to sterile 125 mL screw-capped bottles. Immerse bottles in a 35°C water bath for 5 minutes. Do not directly transfer bottles from the refrigerator to 56°C. Transfer the bottles to a 56°C water bath and begin timing for 30 minutes. To avoid contamination, do not allow the level of the water in the bath to come in contact with the lip of the screw cap. It is best to leave one inch between the serum level in the bottle and the lip of the cap and to fill the water bath to a level just slightly above the level of the serum. To assure even heating of the serum, swirl the bottle(s) every ten minutes. **Note:** Some suppliers provide serum already heat-inactivated.

Transferring the Cell Line

- When the cell line forms a confluent layer, remove all the medium and replace it with 3 mL of Phosphate Buffered Saline (PBS) (ATCC® cat. 30-2200). Incubate T-25 flask at 35°-37°C for 25-30 min.
- Remove all the PBS and replace it with 2 mL of 0.25% (w/v) trypsin dissolved in Hank's Balanced Salt Solution (ATCC® cat. 30-2101).
- Gently distribute the trypsin over the monolayer, remove the trypsin, and place the flask at 35°-37°C for 10 min.
- 4. Add 2 mL of ATCC® 30-2002 with 10% (v/v) HIFBS and detach any cells still adherent by alternately aspirating the medium into a pipette and discharging the contents over the monolayer.
- Distribute the cell suspension in 0.5 mL aliquots to 4 T-25 flasks containing 10 mL fresh ATCC[®] 30-2002 with 10% (v/v) HIFBS.
- 6. Outgas the flask for 10 seconds with a 95% air, 5% CO₂ gas mixture.
- 7. Incubate in a 35-37°C CO₂ incubator with the cap screwed on tightly.

Storage and Culture Initiation of Toxoplasma


Frozen ampules packed in dry ice should either be thawed immediately or stored in liquid nitrogen. If liquid nitrogen storage facilities are not available, frozen ampoules may be stored at or below -70°C for approximately one week. <u>Do not under any circumstance store frozen ampules at refrigerator freezer temperatures (generally -20°C).</u> Storage of frozen material at this temperature will result in the death of the culture.



Product Sheet

Toxoplasma gondii (ATCC[®] 50860™)

Please read this FIRST

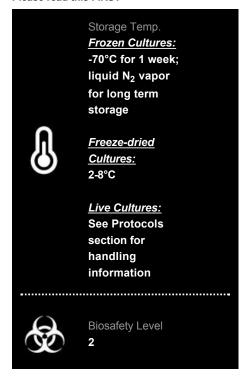

Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: $Toxoplasma\ gondii$ (ATCC® 50860TM)

American Type Culture Collection PO Box 1549 Manassas, VA 20108 USA www.atcc.org

800.638.6597 or 703.365.2700 Fax: 703.365.2750 Email: Tech@atcc.org

Or contact your local distributor

- To thaw a frozen ampule, place it in a 35°C water bath such that the lip of the ampule remains above
 the water line. Thawing time is approximately 2 to 3 minutes. Do not agitate the ampule. Do not leave
 ampule in water bath after it is thawed.
- Immediately after thawing, aseptically transfer contents to a T-25 tissue culture flask containing a fresh monolayer of ATCC® CRL-1634™ cells and 10 mL ATCC® 30-2002 with 3% (v/v) HIFBS.
- 3. Outgas the flask for 10 seconds with a 95% air, 5% CO₂ gas mixture.
- 4. Incubate in a 35-37°C CO₂ incubator with the cap screwed on tightly.

Culture Maintenance

- Remove the medium from a fresh confluent monolayer of CRL-1634™ cells in a T-25 tissue culture flask and replace it with 10 mL of ATCC® 30-2002 with 3% (v/v) HIFBS.
- 2. To transfer the *Toxoplasma* culture, remove the old medium containing the organism and centrifuge at 1300 x g for 10 min.
- Remove all but 0.5 mL of the supernatant and resuspend the cell pellet. Transfer the resuspended pellet to the fresh flask of CRL-1634™ cells.
- 4. Outgas the flask for 10 seconds with a 95% air, 5% CO₂ gas mixture.
- 5. Incubate in a 35°C to 37°C CO₂ incubator with the cap screwed on tightly.



Cryopreservation

Harvest and Preservation

- To harvest the Toxoplasma culture, detach any remaining tissue culture cells (infected and uninfected) by scraping the surface of the flask with a cell scraper.
- 2. Transfer the cell suspension (including parasites) to 15 mL plastic centrifuge tubes. Centrifuge at $1300 \times g$ for 10 min.
- 3. Remove all but 0.5 mL of the supernatant from each tube, resuspend the cell pellets, and pool them to a single tube.
- 4. Pass the resulting cell suspension through a syringe equipped with a 27 gauge 1/2 in needle to break up any remaining cells. Adjust the parasite concentration to 2.0 - 4.0 x 10⁷ cells/mL with fresh medium or PBS. NOTE: If the concentration of parasites is too low, centrifuge at 1300 x g for 10 min and resuspend in the volume of fresh medium or PBS required to yield the desired concentration.
- Prepare a cryoprotective solution containing 15% (v/v) DMSO and 50% (v/v) HIFBS in fresh medium or PBS.
- 6. Mix the cell preparation and cryoprotective solution in equal portions. The final concentration will be 1.0 2.0×10^7 cells/mL, 7.5% DMSO, and 25% HIFBS. The time from the mixing of the cell preparation and cryoprotective solution to the start of the freezing process should be no less than 15 min. and no more than 30 min. NOTE: To prevent culture contamination, penicillin-streptomycin solution (ATCC $^{\textcircled{8}}$ 30-2300) may be added to a final concentration of 50 to 100 I.U./mL penicillin and 50 to 100 µg/mL streptomycin.
- 7. Dispense in 0.5 mL aliquots to 1.0-2.0 mL sterile plastic screw-capped cryules (special plastic vials for cryopreservation).
- 8. Place vials in a controlled rate freezing unit. From room temperature cool at -1°C/min to -40°C. If freezing unit can compensate for the heat of fusion, maintain rate at -1°C/min through heat of fusion. At -40°C plunge ampules into liquid nitrogen. Alternatively, place the vials in a Nalgene 1°C freezing apparatus. Place the apparatus at -80°C for 1.5 to 2 hours and then plunge ampules into liquid nitrogen. (The cooling rate in this apparatus is approximately -1°C/min.)
- 9. Store frozen ampules in either the vapor or liquid phase of a nitrogen refrigerator.
- 10. To thaw a frozen ampule, place it in a 35°C water bath such that the lip of the ampule remains above the water line. Thawing time is approximately 2 to 3 minutes. Do not agitate the ampule. Do not leave ampule in water bath after thawed.
- 11. Immediately after thawing, aseptically transfer contents to a T-25 tissue culture flask containing a fresh monolayer of ATCC® CRL-1634™ cells and 10 mL ATCC® 30-2002 with 3% (v/v) HIFBS.
- 12. Outgas the flask for 10 seconds with a 95% air, 5% CO₂ gas mixture.
- 13. Incubate in a 35-37°C CO₂ incubator with the cap screwed on tightly.



References

References and other information relating to this product are available online at www.atcc.org.



Biosafety Level: 2


Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes



Product Sheet

Toxoplasma gondii (ATCC[®] 50860™)

Please read this FIRST

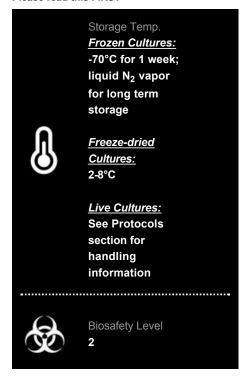

Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: $Toxoplasma\ gondii$ (ATCC® 50860TM)

American Type Culture Collection PO Box 1549 Manassas, VA 20108 USA www.atcc.org

800.638.6597 or 703.365.2700 Fax: 703.365.2750 Email: <u>Tech@atcc.org</u>

Or contact your local distributor

for Health.

ATCC Warranty

The viability of ATCC® products is warranted for 30 days from the date of shipment, and is valid only if the product is stored and cultured according to the information included on this product information sheet. ATCC lists the media formulation that has been found to be effective for this strain. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this strain. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans.

While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, ATCC is not liable for damages arising from the misidentification or misrepresentation of cultures.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this product. The MTA is also available on our Web site at www.atcc.org

Additional information on this culture is available on the ATCC web site at www.atcc.org.

© ATCC 2014. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [08/22]