



Product Sheet

# *Methanolobus oregonensis* (ATCC® BAA-928™)

Please read this **FIRST**



## Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

## Citation of Strain

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Methanolobus oregonensis* (ATCC® BAA-928™)

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor

## Description

**Designation:** OCM 99 [DSM 5435, OGI 99, WAL1]

**Deposited Name:** *Methanolobus oregonensis* (Liu et al.) Boone

## Propagation

### Medium

Methanolophilus oregonese

### Growth Conditions

**Atmosphere:** hydrogen (H<sub>2</sub>), 80%; carbon dioxide (CO<sub>2</sub>), 20%

**Temperature:** 37.0°C

### Propagation Procedure

1. Sterilize the top of the Balch tube by spraying it with 70% ethanol and then flaming the top.
2. If needed, exchange the gas in the test tube for 100% N<sub>2</sub>.
3. If the medium is pink (see discussion about resazurin) add 2.0 ml of reducing agent (3% cysteine stock solution) per 100 ml of medium. Let the medium sit at room temperature for 10 to 20 minutes, until the resazurin becomes colorless, before inoculating.
4. Thaw the frozen vial and using a gassed 1.0 ml syringe tipped with 22-gauge needle, withdraw the cell suspension from the vial and transfer to the Balch tube. Plate 0.1 ml of the inoculated culture onto a non-selective medium and incubate the plate aerobically at 37°C. Incubate broth tube at 37°C.
5. Growth should be detected in the broth within 3 to 4 days. No growth should be detected on the aerobic plate.

### ANAEROBIC CONDITIONS:

- A. Balch tubes (available from Bellco Glass, Vineland, NJ; are specially designed for anaerobic work and use an aluminum crimp cap to hold a rubber stopper in place. Needles can easily be inserted through the stopper, and the tubes can be pressurized to 2 atm. Alternatively, serum vials may be used, or screw cap tubes with butyl rubber stoppers, in the latter case the stopper may be removed and the tube placed under a cannula system that dispenses sterile, oxygen free gas for addition of reducing agents or inoculation.
- B. Resazurin is a commonly used redox indicator that is pink when the redox potential is above 50 mv, and colorless when the redox potential is below 110 mv. i.e. highly reducing. Most strict anaerobes require this low redox potential for optimum growth.
- C. To obtain a fully reduced medium, it is necessary that the medium be anoxic and that a reducing agent be added. Common reducing agents are sodium sulfide, cysteine, dithiothreitol, and titanium citrate.
- D. Syringes can be made anaerobic by one of two methods.

## References

References and other information relating to this product are available online at [www.atcc.org](http://www.atcc.org).

## Biosafety Level: 1

Appropriate safety procedures should always be used with this material. Laboratory safety is discussed in the current publication of the *Biosafety in Microbiological and Biomedical Laboratories* from the U.S. Department of Health and Human Services Centers for Disease Control and Prevention and National Institutes for Health.

## ATCC Warranty

ATCC® products are warranted for 30 days from the date of shipment, and this warranty is valid only if the product is stored and handled according to the information included on this product information sheet. If the ATCC® product is a living cell or microorganism, ATCC lists the media formulation that has been found to be effective for this product. While other, unspecified media may also produce satisfactory results, a change in media or the absence of an additive from the ATCC recommended media may affect recovery, growth and/or function of this product. If an alternative medium formulation is used, the ATCC warranty for viability is no longer valid.

## Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans. While ATCC uses reasonable efforts to include accurate and up-to-date information on this product sheet, ATCC makes no warranties or representations as to its accuracy. Citations from scientific literature and




Product Sheet

## ***Methanolobus oregonensis*** **(ATCC® BAA-928™)**

---

**Please read this FIRST**

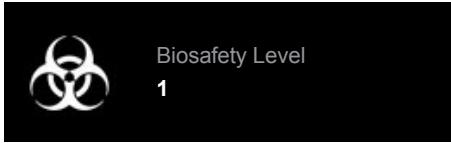

---

**Intended Use**

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

---

**Citation of Strain**

If use of this culture results in a scientific publication, it should be cited in that manuscript in the following manner: *Methanolobus oregonensis* (ATCC® BAA-928™)

patents are provided for informational purposes only. ATCC does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. ATCC is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of materials on deposit, ATCC is not liable for damages arising from the misidentification or misrepresentation of such materials.

Please see the enclosed Material Transfer Agreement (MTA) for further details regarding the use of this product. The MTA is also available on our Web site at [www.atcc.org](http://www.atcc.org)

Additional information on this culture is available on the ATCC web site at [www.atcc.org](http://www.atcc.org).

© ATCC 2013. All rights reserved. ATCC is a registered trademark of the American Type Culture Collection. [02/05]

American Type Culture Collection  
PO Box 1549  
Manassas, VA 20108 USA  
[www.atcc.org](http://www.atcc.org)

800.638.6597 or 703.365.2700  
Fax: 703.365.2750  
Email: [Tech@atcc.org](mailto:Tech@atcc.org)

Or contact your local distributor